




CARFIL QUALITY
YOUR SUPPLIER FOR ALL LABORATORY NEEDS

Ergonomic Table Top Stereo Microscope

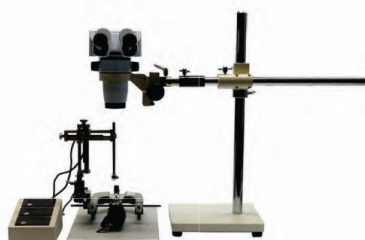

The standard microscope set-up is with the 0,3x objective (ABCD) is designed for use with a stereotaxic instruments because of the 276mm working distance. It can also be used for microsurgery.

The working distance with this unit is 276mm and a magnification upto 15x. If the 0,3x objective is changed for a 0,5 objective (E), the working distance is 126mm and a magnification upto 25x.

Modular Set-up

- Head.
- Zoom Unit.
- Eyepieces.
- Objective.
- Stand.

Specifications

- Ergonomical binocular 0-35° tilting head with 2 eyepieces.
- Large working distance for use during stereotaxic procedures.
- Wide field of view.
- Central Zoom Unit 1:8 .
- 19kg very stable stand with one articulating arm.
- Planachromatic common main 0,3 x objective working distance upto 276mm.
- Optional objective 0,5x for working distance upto 25x.

Ordering number

18002020	Table Top Microscope with 0,3x objective
18004005	Optional 0,5x objective



PAVAN SERVICE BVBA

Beyntellus 3 - B-2360 Oud-Turnhout - Tel. +32 (0)14 45 13 10 - admin@carfil.be - BTW BE 0420.681.872

www.carfil.be